Dear Lumoral user,

Lumoral is intended to benefit your health and well-being by improving your oral hygiene. Use it regularly with Lumorinse to reduce plaque, prevent tooth decay and reduce gingivitis.

We recommend at least 2 times a week for 10 minutes at a time.

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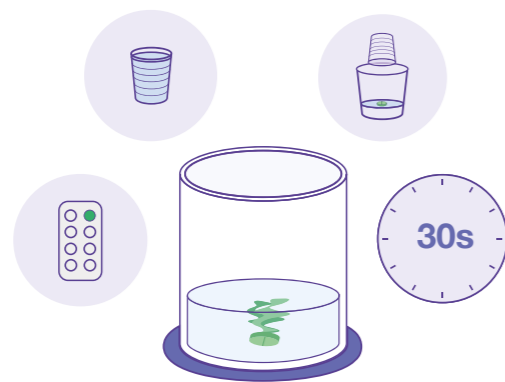
PREPARATIONS


1. Measure **30 ml** of water in a glass or a cup.
2. Take one tablet out of the blister and let the tablet dissolve in the water for **30 s.** or until fully dissolved.

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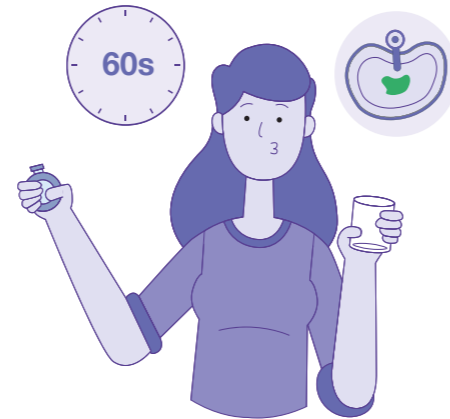SWISH YOUR MOUTH WITH LUMORINSE

1. Swish the mouth rinse for **60 s.**
2. Spit out in the sink.
3. Do not rinse your mouth with water.

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APPLY THE LUMORAL LIGHT

1. Connect the device to the powerbank.
2. Insert the mouthpiece between the lower and upper dental arch, close the mouth gently and press the power button to turn on the light.
3. Keep in mouth until the light goes off (**10 minutes**).

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FINISHING UP


1. Wash the mouthpiece with warm water and let it dry.
2. Brush your teeth with fluoride toothpaste for **2 minutes**. Remember to clean the interdental spaces.
3. Enjoy the feeling of the perfectly clean teeth!

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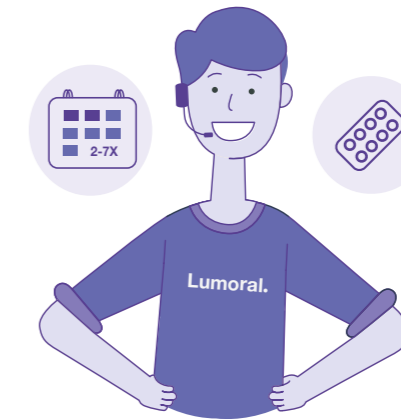USE REGULARLY

1. Repeat every day or at least **2 times** a week.
2. Restock your Lumorinse at **lumoral.com**
3. If you have any questions, we will be happy to help at **info@lumoral.com**

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Koite Health Oy
Kutojantie 2 C
02630 Espoo
Finland
koitehealth.com

Sensations during the use

Lumoral produces heat. You will feel warmth in the teeth and gums during treatment. Heat and near-infrared light may widen the blood vessels (vasodilation). This may create a feeling of weight or pulsing in the gums. Your mouth may temporarily produce more saliva. Saliva protects teeth and has an important oral health-promoting function. Lumorinse does not visibly stain teeth or oral mucosa. If you find stainings, we recommend meticulous brushing of teeth and tongue.

Restrictions of use

Consult your physician, if you have oral candidiasis (oral thrush).

Avoid looking directly at the light source.

Do not connect the mouthpiece to any output, which is powered directly from the electric grid.

Disconnect Lumoral from the powerbank before washing the mouthpiece.

Refer to instruction before use.

Restock your Lumorinse supply at Lumoral.com

Any Questions?
Please check out the FAQ at lumoral.com
Or write to us at info@lumoral.com

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